

**OLYMPUS**<sup>®</sup>

Your Vision, Our Future

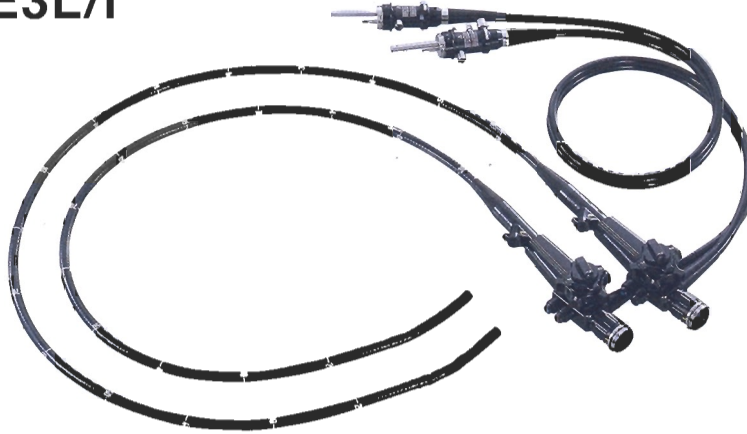
COLONOFIBERSCOPE

**OLYMPUS CF TYPE E3L/I**

**Advanced Fiberoptics Designed for Colonoscopy**

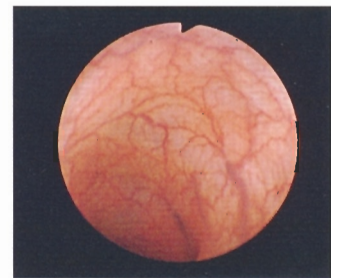
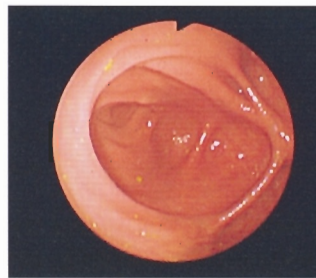


# OLYMPUS CF TYPE E3L/I




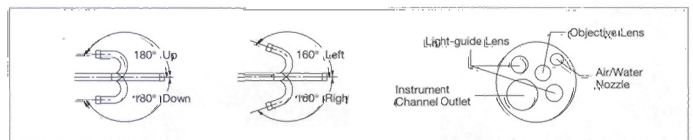
## Expanded 140° field of view and bright image give you clear vision and ease of use

- Wide 140° field of view and extensive 4-way angulation (180° up/down and 160° right/left).
- Both the distal end and insertion tube diameters have been reduced to 13.8 mm and 12.9 mm respectively while keeping the channel diameter the same as its predecessor's.
- Ample 3.2 mm diameter channel comfortably accommodates a wide variety of endotherapy accessories.
- Olympus high-resolution fiberoptics offer good image quality to clearly depict mucosal structures inside the colon.
- All parts of the scope can be easily cleaned and the entire scope can be immersed in a disinfectant solution.
- Ergonomically designed control section facilitates easy handling.



### Specifications

<b>Optical System</b>	Field of view	140°
	Direction of view	0° Forward viewing
	Depth of field	5 to 100 mm
<b>Distal End</b>	Outer diameter	13.8 mm
<b>Insertion Tube</b>	Outer diameter	12.9 mm
<b>Bending Section</b>	Angulation range	Up 180°, Down 180°, Right 160°, Left 160°
<b>Working Length</b>		L : 1680 mm, I : 1330 mm
<b>Total Length</b>		L : 2000 mm, I : 1650 mm
<b>Instrument Channel</b>	Inner diameter	3.2 mm
	Minimum visible distance	2 mm from distal end
	Endotherapy accessory entrance/exit position in field of view	



Specifications, design and accessories are subject to change without any notice or obligation on the part of the manufacturer.

**OLYMPUS**

**OLYMPUS MEDICAL SYSTEMS CORP.**  
Shinjuku Mitsubishi, 3-1 Nishi-Shinjuku 2-chome, Shinjuku-ku, Tokyo 163-0914, Japan  
**OLYMPUS MEDICAL SYSTEMS EUROPA GMBH**  
Postfach 107 49-03, 20034 Hamburg // Wardenstrasse 14-18, 20097 Hamburg, Germany  
**OLYMPUS AMERICA INC.**  
3500 Corporate Parkway, P.O. Box 610, Center Valley, PA 18034-0610, U.S.A.  
**OLYMPUS LATIN AMERICA, INC.**  
5301 Blue Lagoon Drive, Suite 200 Miami, Florida 33126-2097, U.S.A.  
**KEYMED LTD.**  
Keymed House, Stock Road, Southend-on-Sea, Essex SS2 5QH, England

**OLYMPUS SINGAPORE PTE LTD.**  
491B River Valley Road #12-01/04 Valley Point Office Tower, Singapore 248373  
**OLYMPUS HONG KONG AND CHINA LIMITED**  
Room 1520-1527, Ocean Centre, 5 Canton Road, Tsimshatsui, Kowloon, Hong Kong  
**OLYMPUS (BEIJING) SALES & SERVICE CO., LTD.**  
Room 1202, NCII Tower, A12 Jianguomenwai Avenue Chaoyang District Beijing, 100022, China  
**OLYMPUS MOSCOW LIMITED LIABILITY COMPANY**  
117071, Moscow, Malaya Kaluzhskaya 19, Bld. 1, #2, Russia  
**OLYMPUS AUSTRALIA PTY LTD**  
31 Gilby Road, Mount Waverley, VIC, 3149, Australia

www.olympus.com