


OLYMPUS

Your Vision, Our Future

Endo) **Therapy**

SD-990-15/25

Single-use monofilament electro-surgical polypectomy snare.

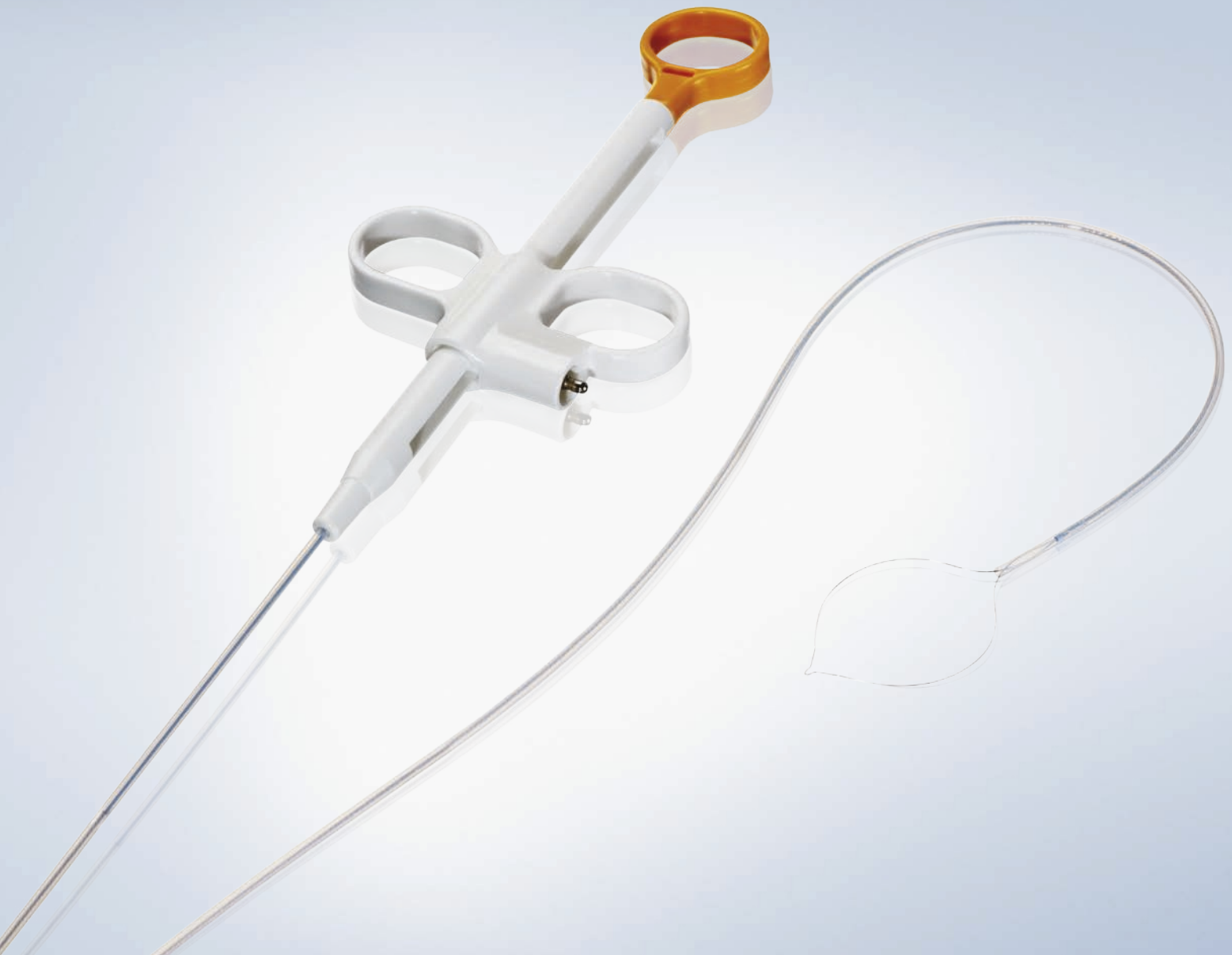

SD-990-15/25

Specifications

SD-990-15, WA95036A	
Loop diameter	15 mm
Snare loop shape	Oval
Working length	2300 mm
Compatible working channel diameter	≥2.8 mm
Outer tube diameter	2.3 mm
Snare wire diameter	0.3 mm
Delivery quantity	10/box

SD-990-25, WA95038A	
Loop diameter	25 mm
Snare loop shape	Oval
Working length	2300 mm
Compatible working channel diameter	≥2.8 mm
Outer tube diameter	2.3 mm
Snare wire diameter	0.3 mm
Delivery quantity	10/box



Features

Hard and rigid wire

The monofilament wire has an outstanding level of stiffness and therefore gives a strong and firm grip. This enables the snare to securely snatch the mucosa without slipping off, facilitating the resection of flat polyps in particular.

Monofilament and thin wire

Due to the 0.3-mm-thin monofilament wire, the snare achieves a very high cutting speed compared to standard snares. The shortened exposure to the tissue means a low coagulation impact on the mucosa and the submucosa.

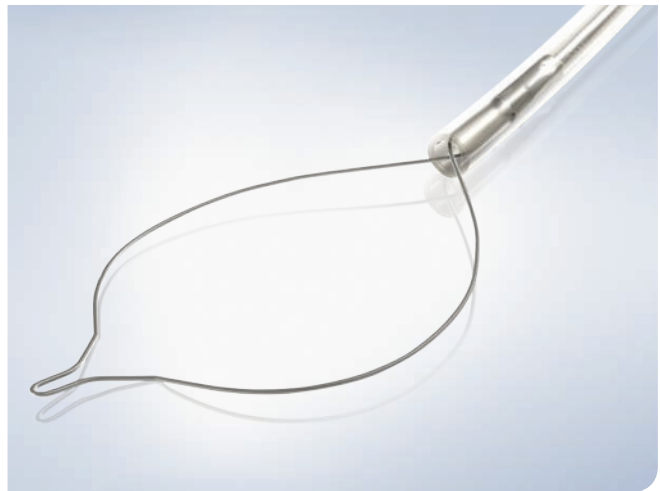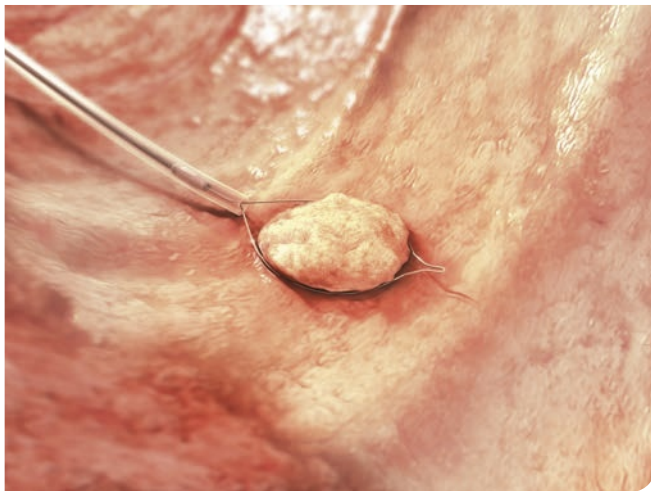

Manufactured by:

Wolfgang Meiners Medizintechnik
Niederstrasse 87 | 40789 Monheim am Rhein | Germany

Distributed by:

Olympus Europa Holding GmbH
Wendenstrasse 14-18 | 20097 Hamburg | Germany

Specifications, design and accessories are subject to change without any notice or obligation on the part of the manufacturer.

OLYMPUS

OLYMPUS EUROPA HOLDING GMBH

Postbox 10 49 08, 20034 Hamburg, Germany
Wendenstrasse 14-18, 20097 Hamburg, Germany
Phone: +49 (0)40 237 730, Fax: +49 (0)40 230 761
www.olympus-europa.com